

# Caas IV-LINQ


*Real time co-registration between angiography and IVUS/OCT*

## Why co-registration?

Co-registration integrates angiography and IVUS/OCT to provide the most detailed view of the lesion related to its exact location in the coronary tree. Co-registration facilitates optimal assessment of the landing zone of the stent.




*Vendor independent: runs with most X-ray, IVUS and OCT systems*

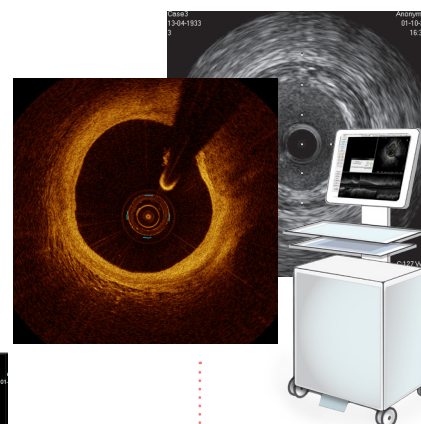

*Only two angiograms needed:*

- one with contrast
- one without contrast

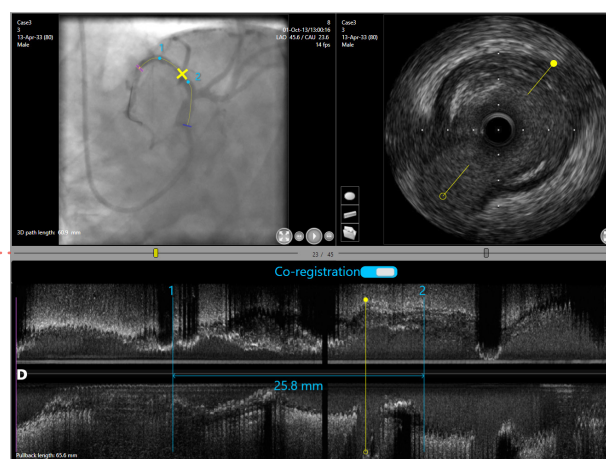

*Using real time streaming, co-registration is ready to be used as soon as the pullback is finished!*

## Key results

- Stent landing zone planned on IVUS/OCT pullback visualized on angiogram.
- Add bookmarks to indicate important vessel sections.
- Diameter and area measurements in cross-sectional view of IVUS and OCT.



Scan QR code for more information or free trial license.

PIE  
MEDICAL  
IMAGING






# Caas IV-LINQ

## Added value of co-registration between angiography and IVUS/OCT

1. Providing exact anatomical registration between angiography and IVUS
2. Co-registration enables non-contrast PCI procedures
3. Precisely identifies locations especially in tortuous vessels
4. Allows visualization of anatomical location and composition of CTO registered to angiogram
5. Facilitates side by side placement of stents

### Added value #4

Co-registration (where the IVUS catheter can cross the occlusion) allows anatomical information of occluded segment on the angiogram.

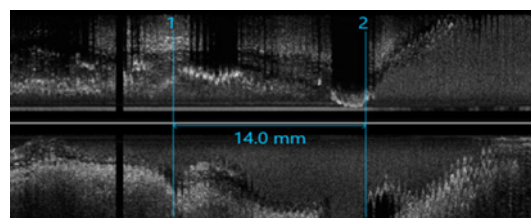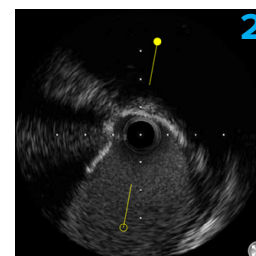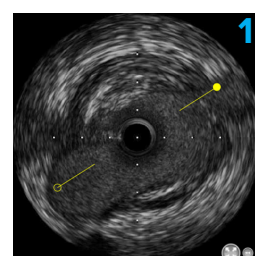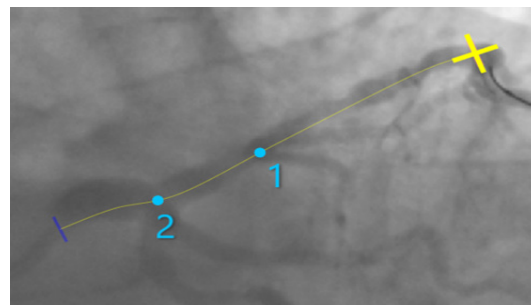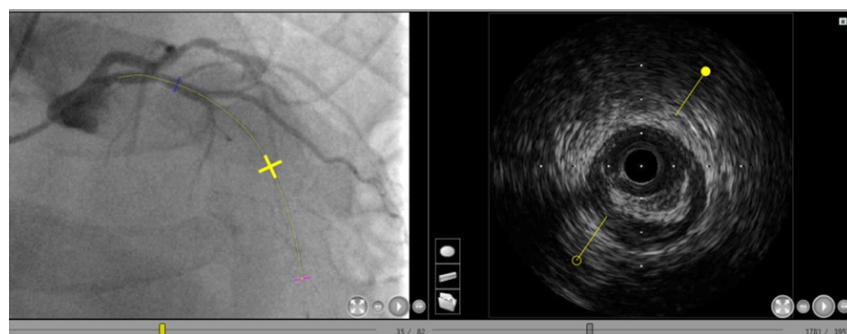

### Added value #1

Angiogram shows two bifurcations at bookmark 1, IVUS shows this as a trifurcation.

**Co-registration appears accurate, robust and induced a minimal delay on the normal cathlab activities.**

*Published in Cardiovascular  
Revascularization Medicine*

### Pie Medical Imaging stands for:

- Having the best understanding of clinical workflows
- Being a knowledge partner where our customers can build on
- Fast and intuitive software solutions

### Quality Assurance:

Pie Medical Imaging develops, produces and sells products in accordance with international accepted standards. CAAS Workstation is FDA 510(k) cleared and MDD compliant (CE marked).

Philipsweg 1  
6227 AJ Maastricht  
The Netherlands

t +31 (0)43 328 13 28  
e [pmi@pie.nl](mailto:pmi@pie.nl)  
i [www.piemedicalimaging.com](http://www.piemedicalimaging.com)

**PIE  
MEDICAL  
IMAGING**